

Redefining Comfort.


**NEW**



# Redefining Comfort.

Comfort is unique. Comfort is personal. Comfort is key to successful CPAP therapy. That's why we believe comfort needs to be measured across a full night's sleep, not just during the initial consultation.

With technologies that work together to provide stability, adaptability and breathability throughout the night, F&P Vitera is a mask that redefines comfort.

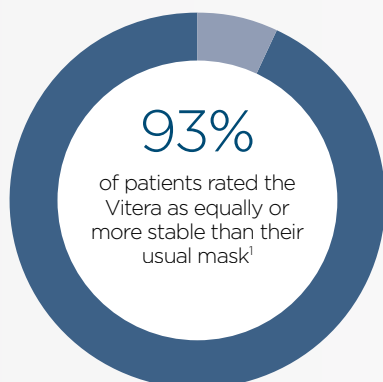

## Stable and secure with dynamic stability.

The stability bar and the RollFit™ XT seal provide dynamic stability to keep the seal comfortable and secure during the night.

## Freedom to move with RollFit XT.


- Dynamically adapts as patients move, helping to keep the seal in place throughout the night
- Soft silicone seal rolls back and forth, minimising pressure on the nose
- Designed to fit a wide range of face and nose shapes
- Seal shape designed to minimise leak into the eyes



88%

of patients could  
sleep in their  
preferred sleep  
position throughout  
the night<sup>2</sup>



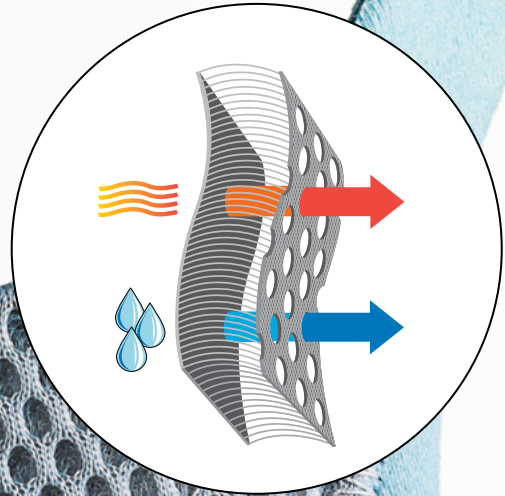


## Keep cool and comfortable with VentiCool.

Sweating during sleep is a well-known symptom of patients with OSA.<sup>4</sup>

Vitera headgear features proprietary breathable fabric that allows more air and moisture to escape than the material used in our leading range of full-face masks.\*

VentiCool technology helps to keep patients cool and comfortable during sleep.

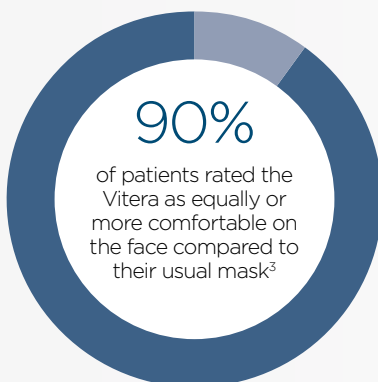

\* Internal testing shows the headgear with VentiCool technology allows 21 times more air-flow and 32% more moisture transfer than the headgear used in our leading full-face masks.<sup>5</sup>


# Key features and benefits

## RollFit XT Seal

Dynamically adapts to minimise leaks and pressure on the bridge of the nose

## VentiCool™ Technology

Breathable headgear designed to allow heat and moisture to dissipate

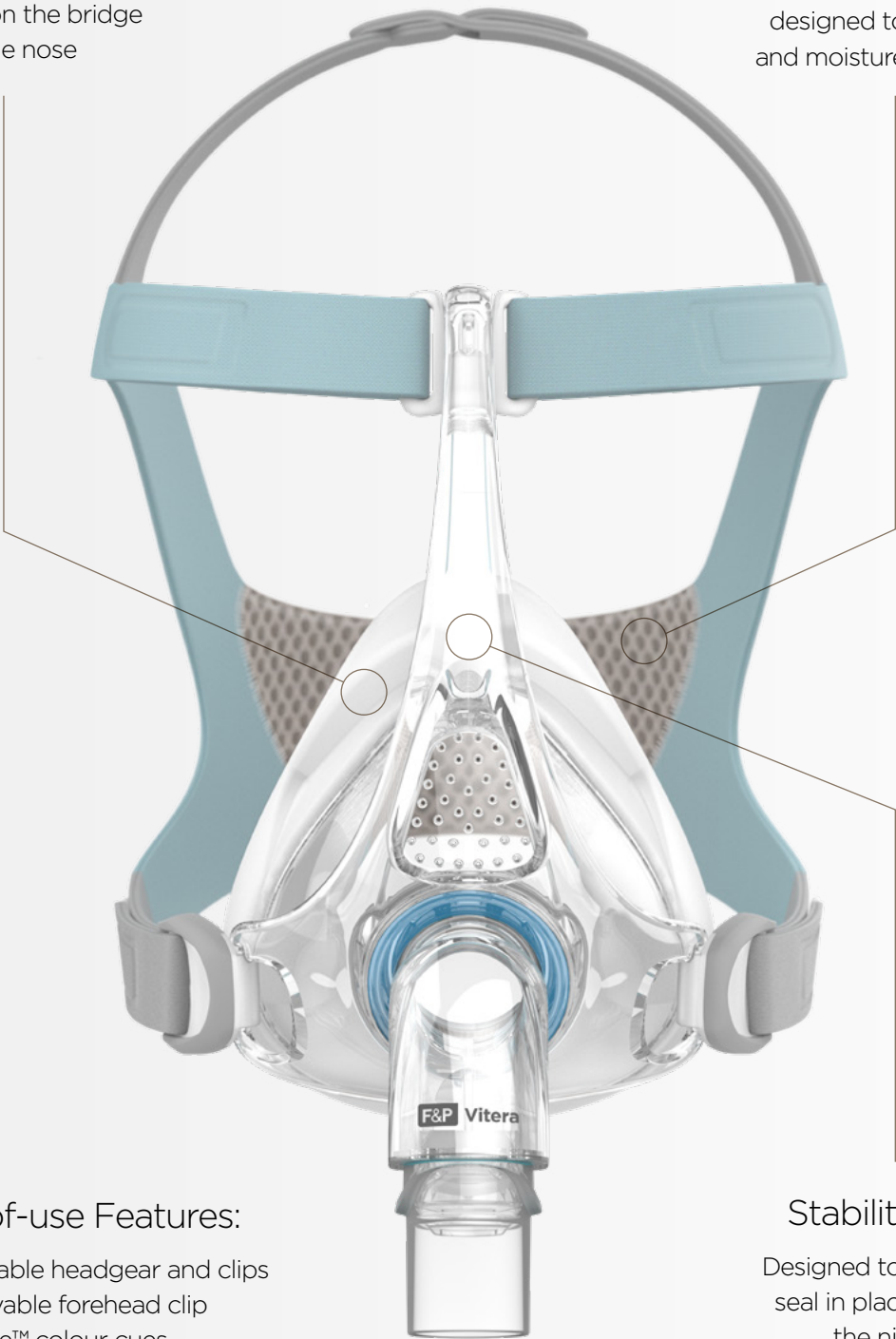

## Ease-of-use Features:

- Adjustable headgear and clips
- Removable forehead clip
- VisiBlue™ colour cues

## Stability Bar

Designed to keep the seal in place during the night

# Fitting guide



**Step 1.**

Hold the front of the mask with one hand and the headgear with the other. Ensure that one headgear clip is unhooked from the frame.



**Step 2.**

Place the seal onto your face and guide the headgear over your head.



**Step 3.**

Hook the unattached headgear clip onto the frame.



**Step 4.**

Gently tighten the headgear straps of the headgear starting with the bottom straps, followed by loosely fitting the blue forehead straps.

Oceania, Canada, Western Europe, Middle East (U Model)

CODE	DESCRIPTION
VIT1SU	Vitera Mask Full Face Small
VIT1MU	Vitera Mask Full Face Medium
VIT1LU	Vitera Mask Full Face Large
VIT1SSU	Vitera Mask Small/Small (two seals)*
VIT1MMU	Vitera Mask Medium/Medium (two seals)*
VIT1LLU	Vitera Mask Large/Large (two seals)*
VIT1MLU	Vitera Mask Medium/Large (two seals)*

Asia, Eastern Europe (Z Model)

CODE	DESCRIPTION
VIT1SZ	Vitera Mask Full Face Small
VIT1MZ	Vitera Mask Full Face Medium
VIT1LZ	Vitera Mask Full Face Large

Spare Parts

CODE	DESCRIPTION
400VIT131	Vitera Mask Full Face Mask No H/G Small*
400VIT132	Vitera Mask Full Face Mask No H/G Medium*
400VIT133	Vitera Mask Full Face Mask No H/G Large*
400VIT111	Vitera Seal Small
400VIT112	Vitera Seal Medium
400VIT113	Vitera Seal Large

CODE	DESCRIPTION
400VIT121	Vitera Headgear Small
400VIT122	Vitera Headgear Medium/Large
400FPH241	F&P Swivel
400VIT151	Vitera Clips and Forehead Clip
400VIT171	Vitera Elbow

\*Not available in all countries.

1. 39 out of 41 participants rated the mask as equally or more stable than their usual mask. Internal validation study conducted on 41 participants in United States of America 2017. 2. 37 out of 41 participants found the mask stable enough so they could sleep in their preferred position throughout the night. Internal validation study conducted on 41 participants in United States of America 2017. 3. 38 out of 41 participants rated the mask as equally or more comfortable on the face compared to their usual mask. Internal validation study conducted on 41 participants in United States of America 2017. 4. Guilleminault, C., and Bassiri, A. Clinical features and evaluation of obstructive sleep apnea-hypopnea syndrome and upper airway resistance syndrome. In: Kryger, M.H., Roth, T., and Dement, W.C. (Eds). Principles and Practice of Sleep Medicine. Elsevier Sanders, Philadelphia, PA, 2015; 1043-1052. 5. Internal study of F&P Vitera testing breathability of fabric with VentiCool technology against F&P Simplus™ and F&P FlexiFit™ headgear materials.

Vitera, RollFit, VentiCool, Simplus, FlexiFit and VisiBlue are trademarks of Fisher & Paykel Healthcare Limited. For patent information, see [www.fphcare.com/ip](http://www.fphcare.com/ip)