
DE Gebrauchsanweisung **EN** Instructions for Use **FR** Mode d'emploi
NL Gebruiksaanwijzing **IT** Istruzioni d'uso **TR** Kullanma Kılavuzu
ES Manual de instrucciones

WM 65860c 11/2020 DE, EN, FR, NL, IT, TR, ES (MX)

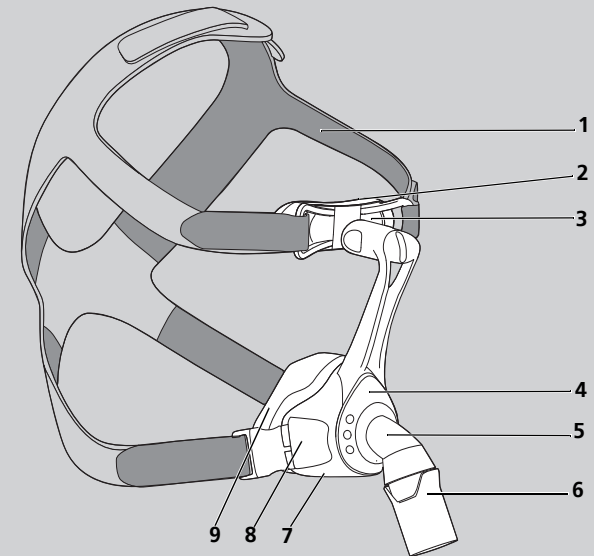
CE 0197

**Löwenstein Medical
Technology GmbH + Co. KG**
Kronsaalsweg 40
22525 Hamburg, Germany
T: +49 40 54702-0
F: +49 40 54702-461
www.loewensteinmedical.com

JOYCEeasy







Nasal Mask
vented, non vented

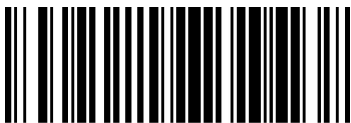WM 65860c

LÖWENSTEIN
medical

LÖWENSTEIN
medical

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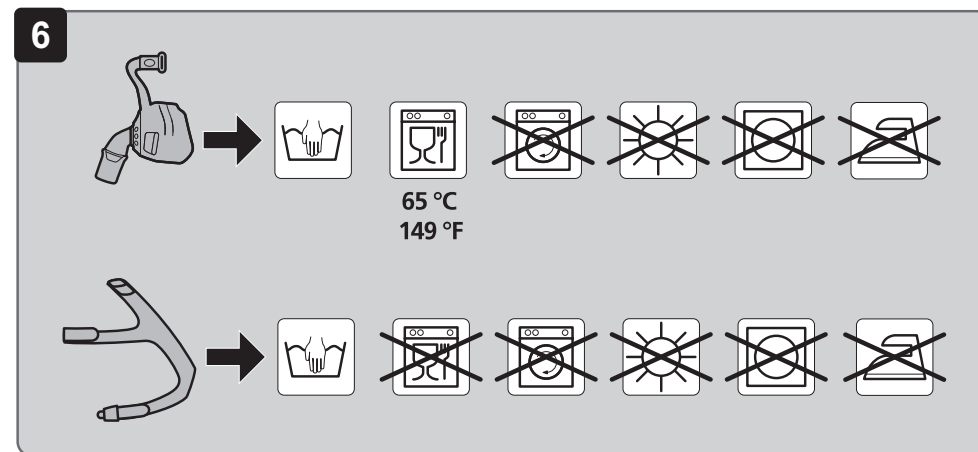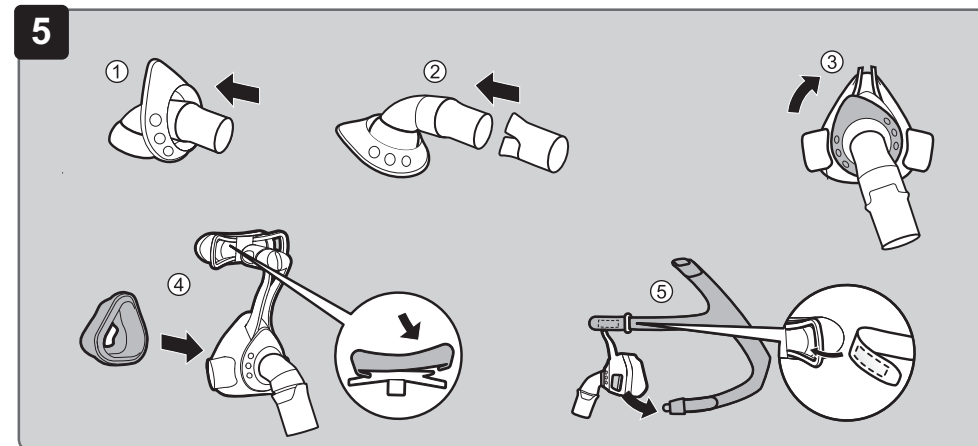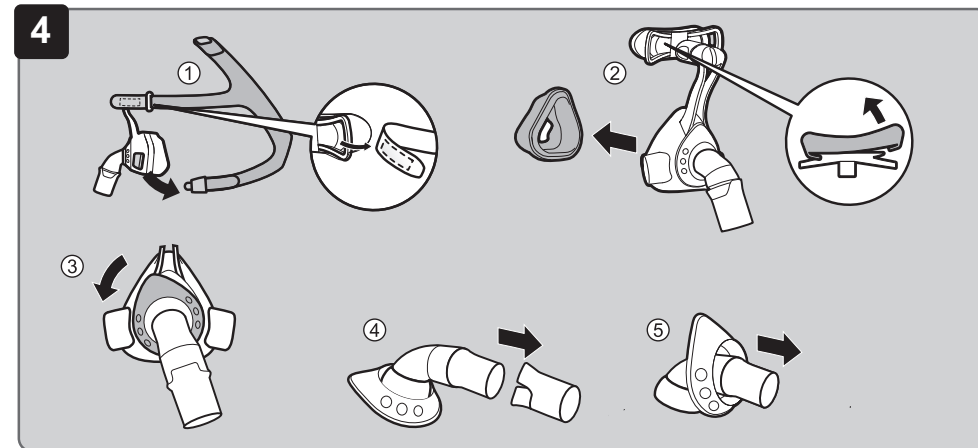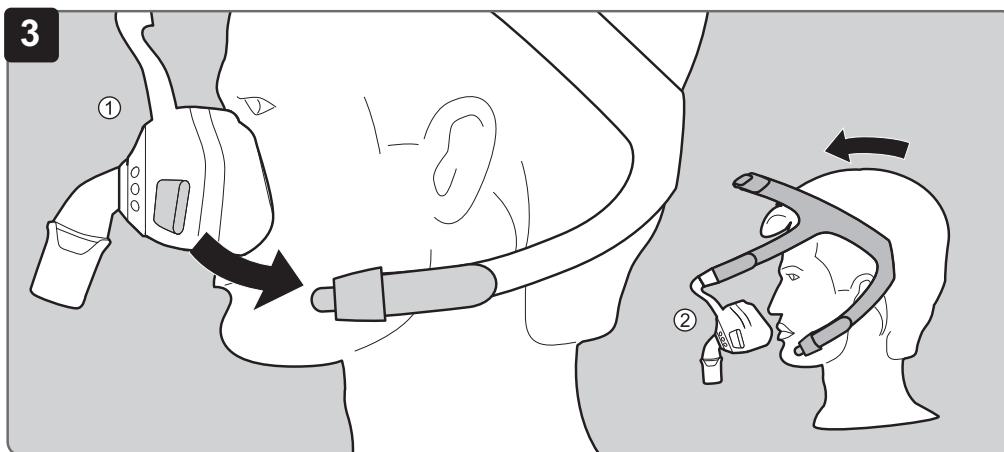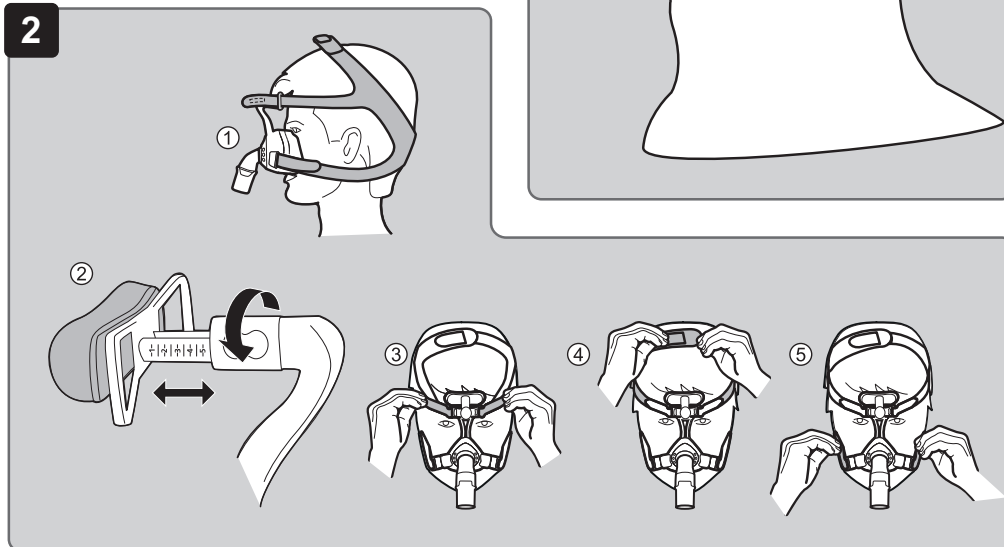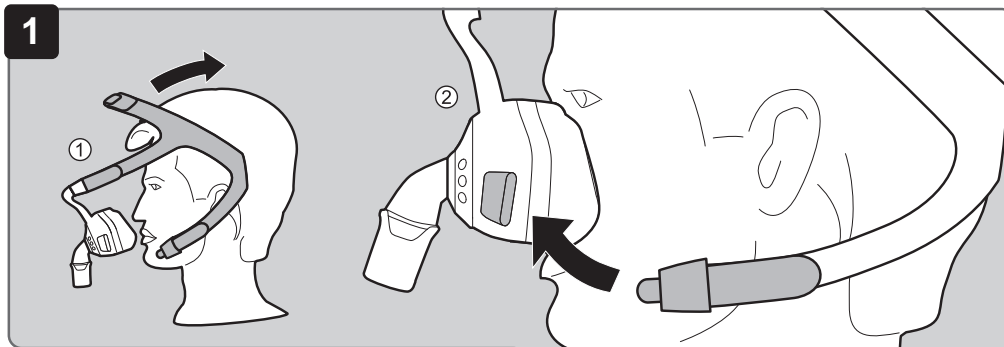1 Operation

The following figures show you how to adjust, apply, remove, dismantle and assemble the mask:

- 1 Putting on the mask
- 2 Adjusting the mask
- 3 Removing the mask
- 4 Dismantling the mask
- 5 Assembling the mask

2 Introduction

2.1 Intended use

The JOYCEeasy mask is used for treating sleep apnea and for non-invasive and non-life-sustaining ventilation of patients with respiratory insufficiency. It serves as a connecting element between the patient and the therapy device.

2.2 Contraindications

The mask may not be used, or may be used only with particular caution, if any of the following symptoms are present:

erosions and ulcerations, skin allergies, flush of the facial skin, pressure marks on the face, claustrophobia, anxiety, deformities of the face or nasopharynx, intake of medication which may induce vomiting, the necessity for immediate intubation. Please also observe the contraindications in the instructions for use of your therapy device.

2.3 Side effects

Nasal congestion, dry nose, dry mouth in the morning, feeling of pressure in the sinuses, irritated conjunctiva, skin rashes, pressure marks on the face, irritating noises when breathing.

3 Safety

Risk of injury from the supply of oxygen!

Oxygen can become deposited in clothing, bedlinen and hair. In conjunction with smoking, naked flames or electrical equipment, it can cause fires and explosions.

- ⇒ Do not smoke.
- ⇒ Avoid naked flames.

Risk of injury from re-inhaling CO₂!

If the mask is used incorrectly, CO₂ may be re-inhaled.

- ⇒ Do not close off the exhalation system of the mask.
- ⇒ Only put on the mask for an extended period if the therapy device is running.

- ⇒ Only use the mask within the quoted therapy pressure range.
- ⇒ Do not use the mask on patients who are unable to take it off themselves.

Risk of injury if the mask slips!

If the mask slips or falls off, the therapy is ineffective.

- ⇒ Monitor patients with restricted spontaneous respiration.
- ⇒ Activate low pressure/leakage alarms on the therapy device.

Risk of injury from anesthetic gases!

Anesthetic gas may escape through the exhalation valve and put third parties at risk.

- ⇒ Never use the mask during anesthesia.

Risk of injury from lack of cleaning!

The mask may show contamination which may put the patient at risk.

- ⇒ Clean the mask before using for the first time (see section entitled "Cleaning and hygiene treatment").
- ⇒ Clean mask regularly.

4 Product description

4.1 Overview

A diagram of the individual parts can be found on the title page.

1. Headgear
2. Forehead cushion
3. Forehead support
4. Retaining ring
5. Elbow
6. Rotating sleeve
7. Mask body
8. Headgear clip
9. Mask cushion

4.2 Compatible devices

With many combinations of device, actual pressure in the mask will not correspond to prescribed therapy pressure, even if the therapy device is displaying the correct pressure. Have the combination of devices adjusted by a physician or specialist dealer so that actual pressure in the mask corresponds to therapy pressure.

4.3 Exhalation system

The mask has an integrated exhalation system. The retaining ring and mask body are shaped so that

there is a gap between these parts. The exhaled air can escape through this gap.

5 Cleaning and hygiene treatment

WARNING

Risk of injury from inadequate cleaning!
Residues may congest the mask, impair the integrated exhalation system and jeopardize therapy success.
⇒ For patients with a compromised immune system or particular background of illness, disinfect mask components daily following consultation with the physician.

5.1 Clean mask

- 1. Dismantle mask (see Figure 4).
- 2. Clean mask in accordance with the table below:

Action	Daily	Weekly
Wash mask components with warm water and a mild detergent.	X	
When washing mask parts, clean thoroughly with a cloth or a soft brush. Or: Put the mask parts in the dishwasher.		X
Wash headgear by hand.		X

- 3. Rinse all parts with clear water.
- 4. Allow all parts to air-dry.
- 5. Perform a visual inspection.
- 6. If necessary: replace damaged parts.
- 7. Re-assemble mask (see Figure 5).

Discolorations of mask components do not impair the functionality of the mask.

5.2 Hygiene treatment (clinical sphere)

In the clinical sphere only: You must subject the mask to a hygiene treatment in the event of a change of patient. Information on the hygiene treatment with a change of patient can be found in a brochure available on the manufacturer's website. We will also send you this brochure on request.

6 Disposal

You can dispose of all parts in domestic waste.

7 Troubleshooting

Fault	Cause	Remedy
Pain from pressure on the face.	Mask too tight.	Loosen headgear slightly.
Draft in the eye.	Mask too loose.	Tighten headgear slightly.
	Mask does not fit.	Contact your specialist dealer.
Therapy pressure is not reached.	Mask not correctly adjusted.	Re-adjust mask (see Figure 2).
	Mask cushion amaged.	Replace mask cushion.
	Patient circuit leaking.	Check connectors, check tubes properly located.

8 Technical specifications


Product class to directive 93/42/EEC	IIa
Dimensions (W x H x D) Size S Size M Size L	100 mm x 100 mm x 80 mm 105 mm x 110 mm x 80 mm 110 mm x 120 mm x 80 mm
Weight Size S Size M Size L	94 g 97 g 98 g
Therapy pressure range	4 hPa - 30 hPa
Tube connection: tapered connection to EN ISO 5356-1	Ø 22 mm (male)
Temperature range: Operation Transport and storage	+5 °C to + 40 °C -20 °C to +70 °C
Flow resistance at 50 l/min at 100 l/min	0.3 hPa 1.0 hPa
Quoted two-figure noise emission value to ISO 4871: - Sound pressure level - Sound power level - Uncertainty factor	16 dB(A) 24 dB(A) 3 dB(A)
Service life	5 years
Service life	up to 12 months ¹
Standards applied	EN ISO 17510-2:2009

¹ The mask materials deteriorate if exposed to e.g. aggressive detergents. In individual cases it may be necessary to replace mask parts sooner.

The right to make design modifications is reserved.

9 Pressure/flow curve

The characteristic pressure/flow curve shows the leakage flow as a function of therapy pressure.

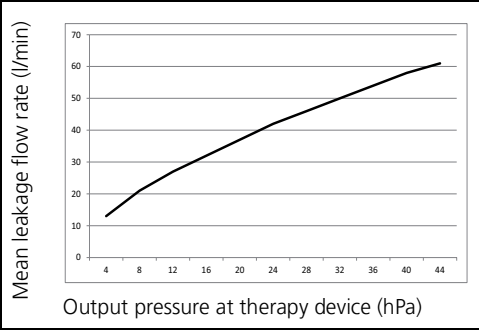10 Materials

Only use the mask after consulting your physician if you are allergic to any of these substances:

Mask cushion, forehead cushion	SI (Silicone)
Retaining ring (hard component)	PP (polypropylene)
Retaining ring (soft component)	TPE (thermoplastic elastomer)
Headgear clip	PA (polyamide), POM (polyoxymethylene)
Forehead support, rotating sleeve, mask body, elbow	PA (polyamide), PC (polycarbonate)
Headgear	EL (Elastane), PES (polyester), PU (polyurethane), UBL loop, CO (cotton), PA (polyamide)

No parts of the mask contain latex, PVC (polyvinyl chloride) or DEHP (diethylhexyl phthalate).

11 Warranty

Löwenstein Medical Technology gives the customer a limited manufacturer warranty on new original Löwenstein Medical Technology products and any replacement part fitted by Löwenstein Medical Technology in accordance with the warranty conditions applicable to the product in question and in accordance with the warranty periods from date of purchase as listed below. The warranty conditions are

available on the website of the manufacturer. We will also send you the warranty conditions on request.

Please bear in mind that any claim to warranty and liability shall be void if neither the accessories recommended in the instructions for use nor genuine replacement parts are used.

In the event of a claim under warranty, contact your specialist dealer.

Product	Warranty periods
Masks including accessories	6 months

12 Declaration of conformity

The manufacturer Löwenstein Medical Technology GmbH + Co. KG (Kronsaalsweg 40, 22525 Hamburg, Germany) hereby declares that the product complies with the relevant regulations of Directive 93/42/EEC governing medical devices. The unabridged text of the Declaration of Conformity can be found on the manufacturer's website.

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