


FlexiFit™ 406

PETITE NASAL MASK

FlexiFit™ 407

FlexiFit™ 406

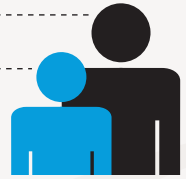

Specification Sheet

A Nasal Delivery Mask




Specifically designed to fit finer featured faces. Incorporating FlexiFit™ Technology with auto-contouring properties, the FlexiFit™ 406 Nasal Mask comes with a unique Glider™ Strap mechanism for maximum comfort and stability. Other unique features ensure a premium sleep experience.

SleepStyle™
SLEEP PERFORMANCE FOR
AN ENERGIZED LIFESTYLE

Fisher & Paykel
HEALTHCARE

Nasal Delivery - *FlexiFit™ 406*


FlexiFit™ 407
FlexiFit™ 406

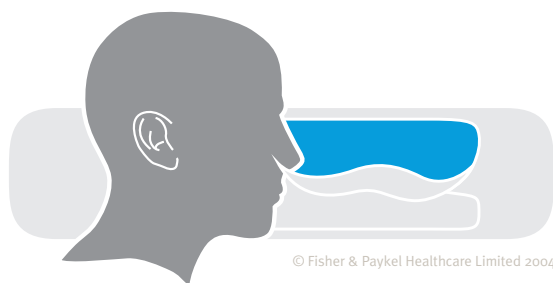

© Fisher & Paykel Healthcare Limited 2004


The FlexiFit™ 406 Petite Nasal Mask offers an over-the-nose mask solution specifically designed to fit finer featured faces. The Petite design is small, light and easy to use. The Stretchgear™ Headgear provides even greater comfort and a maintenance-free bias diffuser ensures quieter operation. The FlexiFit™ 406 Petite Nasal Mask comes as an Out-of-Box Solution fully assembled and ready for use.

UNIQUE MASK FEATURES FOR GREATER FREEDOM, COMFORT AND COMPLIANCE

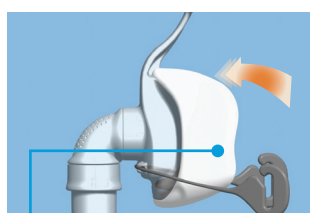


FlexiFit™ Technology Auto-contouring fit

With Foam Cushion
FlexiFit™ Technology automatically adjusts to a wide range of nasal and facial contours with simple one-step ease of fitting and no complicated adjustments.

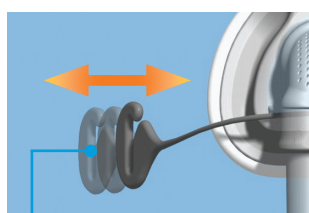

Glider™ Strap Anti-leak system

The sliding Glider™ Strap provides freedom of movement while maintaining a seal. Mask force remains equal when the head rotates from side to side, minimizing the occurrence of leaks.

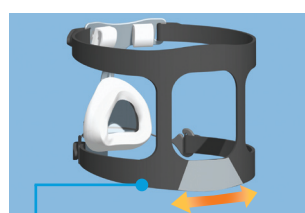


Stretchgear™ Headgear Deluxe mask suspension system

This is a special elastic headgear segment that maintains an equal mask force when the head tilts up or down.




Premium Frosted Silicone Seal Optimizing seal technology

This seal technology has a soft, smooth feel to slide onto the face, and its frosted texture optimizes the mask seal.

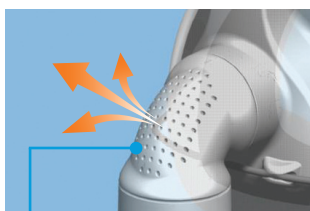


Efficient Bias Diffuser Quiet usage

Maintenance Free
This advanced air diffuser system prevents air from blowing onto the sleeping partner and decreases noise that could potentially disrupt sleep.

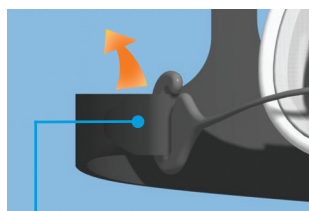

Glider™ Clip One-step removal

Simple one-step removal of headgear. No ongoing headgear adjustments required.

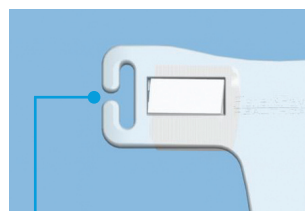

Hidden Strap Release One-step removal

Easy one-step removal of headgear. Enables reassembly without having to make fitting adjustments.

CONTENTS

- Complete mask with headgear assembled and ready for use
- Single one size Silicone Seal designed to fit most finer featured faces
- User friendly Fitting Guide

REORDER CODE: HC406



FOR MORE INFORMATION, PLEASE CONTACT YOUR LOCAL FISHER & PAYKEL HEALTHCARE REPRESENTATIVE:

Manufacturer Fisher & Paykel Healthcare Ltd, 15 Maurice Paykel Place, East Tamaki, Auckland 2013, PO Box 14 348, Panmure, Auckland 1741, New Zealand
Tel: +64 9 574 0100 Fax: +64 9 574 0158 Email: info@fphcare.com Web: www.fphcare.com

Australia Tel: +61 3 9879 5022 Fax: +61 3 9879 5232 **Austria** Tel: 0800 29 31 23 Fax: 0800 29 31 22 **Benelux** Tel: +31 40 216 3555 Fax: +31 40 216 3554 **China** Tel: +86 20 3205 3486 Fax: +86 20 3205 2132 **France** Tel: +33 1 6446 5201 Fax: +33 1 6446 5221 **Germany** Tel: +49 7182 93777 0 Fax: +49 7182 93777 99 **India** Tel: +91 80 4123 6041 Fax: +91 80 4123 6044 **Irish Republic** Tel: 1800 409 011 **Italy** Tel: +39 06 7839 2939 Fax: +39 06 7814 7709 **Japan** Tel: +81 3 3661 7205 Fax: +81 3 3661 7206 **Middle East** Tel: +9 71 4 360 0849 Fax: +9 71 4 368 8110 **Northern Ireland** Tel: 0800 132 189 **Spain** Tel: +34 902 013 346 Fax: +34 902 013 379 **Sweden** Tel: +46 8 564 76 680 Fax: +46 8 36 63 10 **Switzerland** Tel: 0800 83 47 63 Fax: 0800 83 47 54 **UK (EU Authorised Representative)** Fisher & Paykel Healthcare UK Ltd, Unit 16, Cordwallis Park, Clivemont Road, Maidenhead, Berkshire SL6 7BU, UK Tel: +44 1628 626 136 Fax: +44 1628 626 146 **USA** Tel: 1 800 446 3908 or +1 949 453 4000 Fax: +1 949 453 4001

© Fisher & Paykel Healthcare Ltd